


INSTRUCTIONS FOR USE

Dupilumab (200 mg) single-use solution for injection in a Pre-filled Syringe with Needle Shield (PFS-S)

The parts of the Dupixent Syringe are shown in this picture.

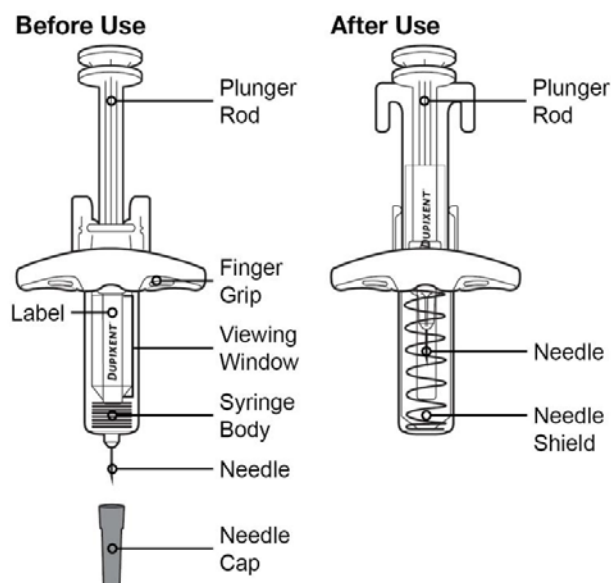

Important information

This device is a single-use Pre-filled Syringe (called “Syringe” in these instructions). It contains 200 mg of Dupixent for injection under the skin (subcutaneous injection).

It is important that you do not try to give yourself or someone else the injection unless you have received training from your healthcare provider.

- To reduce the risk of accidental needle sticks, each Pre-filled Syringe has a needle shield that is automatically activated to cover the needle after you have given your injection.

Do

- ✓ Keep these instructions for future use. Ask your healthcare provider or call sanofi at 1800 818 806 (Australia), if you have further questions.

- ✓ Read all of the instructions carefully before using the Syringe.
- ✓ Ask your healthcare provider how often you will need to inject the medicine.
- ✓ Ask your healthcare provider to show you the right way to use the Syringe before you inject for the first time.
- ✓ Rotate the injection site each time you inject.
- ✓ Keep the Syringe(s) out of the reach of children.
- ✓ Keep unused Syringes in the original carton and store in the refrigerator between 2°C and 8°C.
- ✓ Remove the Syringe from the refrigerator at least 30 minutes before your injection so that it reaches room temperature.


Do not

- ✗ Do not use the Syringe if it has been dropped on a hard surface or damaged.
- ✗ Do not use the Syringe if the Needle Cap is missing or not securely attached.
- ✗ Do not touch the Plunger Rod until you are ready to inject.
- ✗ Do not inject through clothes.
- ✗ Do not get rid of any air bubbles in the Syringe.
- ✗ Do not pull back on the Plunger Rod at any time.
- ✗ Do not re-use the Syringe.
- ✗ Do not keep Dupixent at room temperature for more than 14 days.
- ✗ Do not shake the Syringe at any time.
- ✗ Do not heat the Syringe.
- ✗ Do not freeze the Syringe.
- ✗ Do not put the Syringe into direct sunlight.

STEP 1: REMOVE

Remove the Syringe from the carton by holding the middle of the Syringe Body.

- **Do not pull off the Needle Cap until you are ready to inject.**
- **Do not use the Syringe if it has been dropped on a hard surface or damaged.**

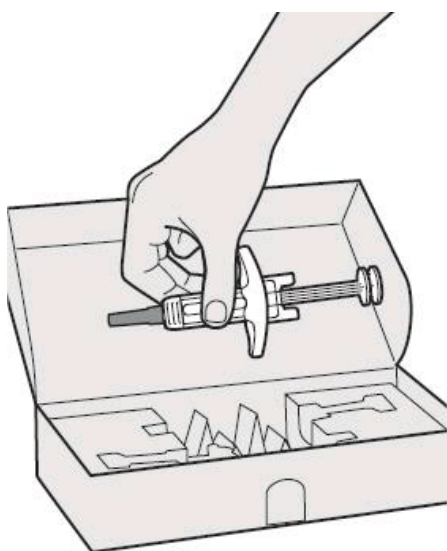

STEP 2: PREPARE


Ensure you have the following:

- the Dupixent Pre-filled Syringe
- 1 alcohol wipe*
- 1 cotton ball or gauze*
- A sharps container* (see Step 12)

* Items not included in the carton

Look at the label:


- Check the expiration date.
- Check that you have the correct product and dose.
- **Do not use the Syringe if the expiration date has passed.**

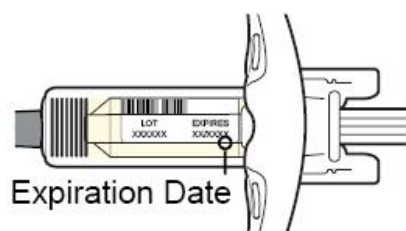

STEP 3: INSPECT

Look at the medicine through the Viewing Window on the Syringe:


- Check if the liquid is clear and colorless to pale yellow.
- Note: You may see an air bubble; this is normal.
- **Do not use the Syringe if the liquid is discolored or cloudy, or if it contains visible flakes or particle.**

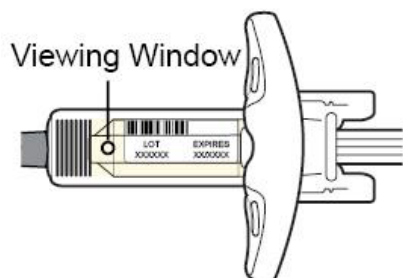

STEP 4: WAIT 30 MINUTES

Lay the Syringe on a flat surface and let it naturally warm to room temperature for at least 30 minutes.


- **Do not heat the Syringe.**
- **Do not put the Syringe into direct sunlight.**
- **Do not keep Dupixent at room temperature for more than 14 days.**

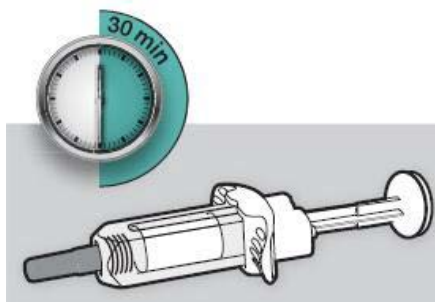

STEP 5: CHOOSE

Select the injection site.

- You can inject into your thigh or stomach, except for the 5 cm (2 inches) around your navel.
- If somebody else gives you the injection, you can also use the upper arm.
- **Do not inject into skin that is tender, damaged or has bruises or scars.**

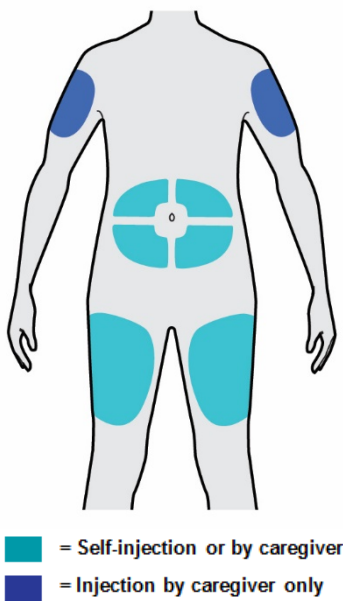


STEP 6: CLEAN

Wash your hands.

Clean the injection site with an alcohol wipe.

Let your skin dry before injecting.

Do not touch the injection site again or blow on it before the injection.

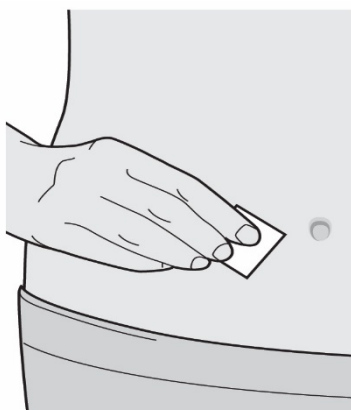


STEP 7: PULL

Hold the Syringe in the middle of the Syringe Body with the Needle pointing away from you and pull off the Needle Cap.


- **Do not put the Needle Cap back on.**
- **Do not touch the Needle.**

Inject your medicine immediately after removing the Needle Cap.

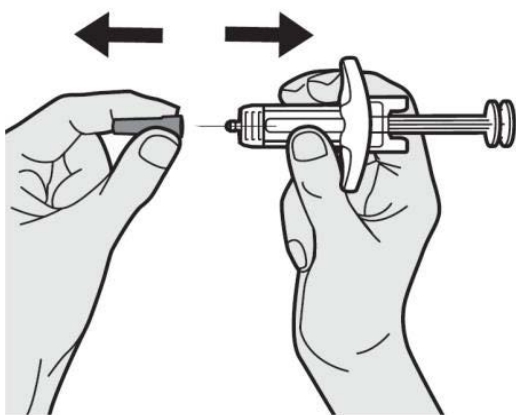


STEP 8: PINCH

Pinch a fold of skin at the injection site, as shown in the picture.

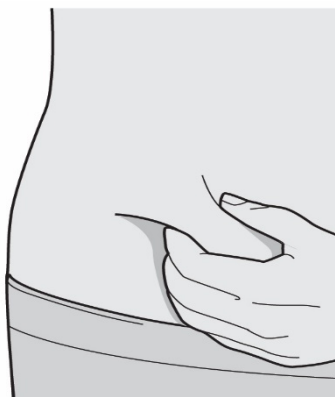

STEP 9: INSERT

Insert the Needle into the fold of skin at roughly a 45° angle.

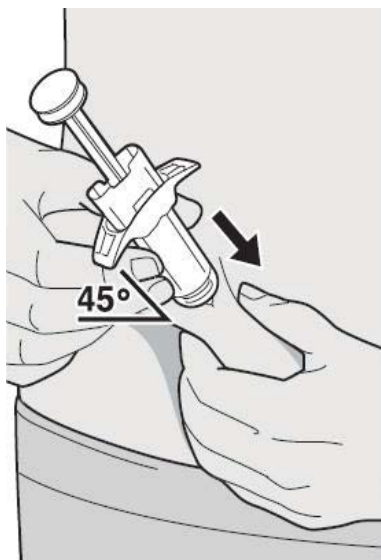

STEP 10: PUSH

Relax the pinch.

Push the Plunger Rod down slowly and steadily as far as it will go until the Syringe is empty.

- Note: You will feel some resistance. This is normal.




STEP 11: RELEASE AND REMOVE

Lift your thumb to release the Plunger Rod until the Needle is covered by the Needle Shield and then remove the Syringe from the injection site.

Lightly press a cotton ball or gauze on the injection site if you see any blood.

- **Do not put the Needle Cap back on.**
- **Do not rub your skin after the injection.**

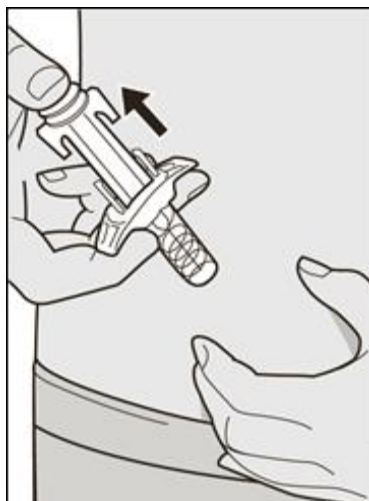

STEP 12: DISPOSE

Dispose of the Syringe and the Needle Cap in a sharps container.

Always keep the container out of the reach of children.

- **Do not put the Needle Cap back on.**
- **Do not throw away the used Syringe in the household rubbish.**

